

WanSen X-RAY TUBE**XD5-0.5-100 (12)****XD5-0.5-100S (12)****XD5-0.5-100SB (12)****XD5-0.5-100SBR(12)****Stationary Anode X-Ray Tube**

- ◆ Designed for DC applications with high input power

This X-ray tube has a 0.5focal point.

The maximum tube voltage that can be used is 100kV.

- ◆ Mounted in the same housing as the high voltage transformer.
- ◆ WanSen product versions comply with IEC standards.

General data

Electrical:

Circuit:

High Voltage Generator constant high voltage generator

Grounding Methods Center-grounded

Nominal X-ray Tube Voltage 100kV

Nominal Value Of Focus :

Focus 0.5 (IEC336)

Nominal Anode Input Power (at 0.1 seconds):

Focus 1.6kW

Nominal Photographic Anode Input Power (at 0.1 seconds):

Focus 1.5kW

Exposure Duty Cycle 1:60 or More

(Exposure Time : Interval Time)

Continuous Anode Input Power 200W

Air Specific Release Kinetic Energy Rate Not Less Than 0.4mGy/s
(75kV, 20mA, 1s, 49cm from focal point (16cm from insulating oil))

Mechanical:

Dimensions:

Overall Length	see dimensional outline
XD5-0.5-100 (12)	135mm
XD5-0.5-100S (12)	143mm
XD5-0.5-100SB (12)	143mm
XD5-0.5-100SBR (12)	155mm
Maximum Diameter	see dimensional outline
XD5-0.5-100 (12)	46mm
XD5-0.5-100S (12)	54mm
XD5-0.5-100SB (12)	58mm
XD5-0.5-100SBR (12)	58mm

Anode Target:

Anode Angle	5°
Target Surface Materials	Tungsten
Inherent Filtration	At Least 0.8mmAl equivalent at 75kV
X-ray Coverage	189×189mm at SID 550mm

Approximate weight:

XD5-0.5-100 (12)	0.5kg
XD5-0.5-100S (12)	0.58kg
XD5-0.5-100SB (12)	0.88kg
XD5-0.5-100SBR (12)	1.0kg

Cooling methodOil immersed (60°C Max.) and convection oil cooling

Tube Holding:

XD5-0.5-100 (12)	Screws for fixing the anode end and cathode end of the glass housing or anode rod
XD5-0.5-100S (12) /XD5-0.5-100SB (12) /XD5-0.5-100SBR (12)	Holding the insulation cylinder

Absolute maximum and minimum ratings

(At any time, these values must not be exceeded.)

Maximum Tube Voltage.....100kV

Minimum Tube Voltage 50 kV

Maximum Tube Current32mA

Maximum Filament Current

Focus3.0A

Filament Voltage (At maximum filament current)

Focus4.5± 0.5V

Thermal Characteristics:

Anode Heat Storage Capacity20KJ (27KHU)

Maximum Anode Heat Dissipation Rate210W

Environmental Limits

Operating Limits (in dielectric oil)

Oil Temperature 10 ~ 60℃

Oil pressure70 ~ 106 kPa

Shipping and Storage Limits:

Temperature -40 ~ 70℃

Humidity10 ~ 90 %

(No condensation)

Atmospheric Pressure50 ~ 106 kPa


Cautions

Read through this page before using the tube.


Since X-ray tube will emit X-rays when it is energized with high voltage, special knowledge is required to handle it. The items below show general cautions for the tube handling.

1. The tube shall be handled or operated only by qualified personnel.
Only a specialist with knowledge of X-ray tube should assemble, maintain and remove the tube.
2. The tube envelope is made of glass. In transporting and handling, sufficient care should be taken not to give strong impact or vibration to the tube.
3. Radiation protection of the tube unit assembled with this tube must be sufficiently taken. And the leakage technique factor of the tube unit must not exceed maximum anode cooling rate of this tube.
4. Regulations and standards require the minimum source-skin distance (SSD) and the minimum filtration of the useful beam.
Use the tube after fulfilling them.
5. The tube might be broken due to only one overload operation.
Provide proper overload protection circuit.
Operate the tube by selecting a proper input condition according to the conditions for operation and tube characteristics charts.
6. When any abnormalities are found in using this tube, immediately switch off the power supply and contact WanSen Service sector.
7. The charts of this specification are indicating standard values.
For usage not described here or for any unclear items, contact WanSen Service sector.


Maximum Rating Charts
(Absolute maximum rating charts)

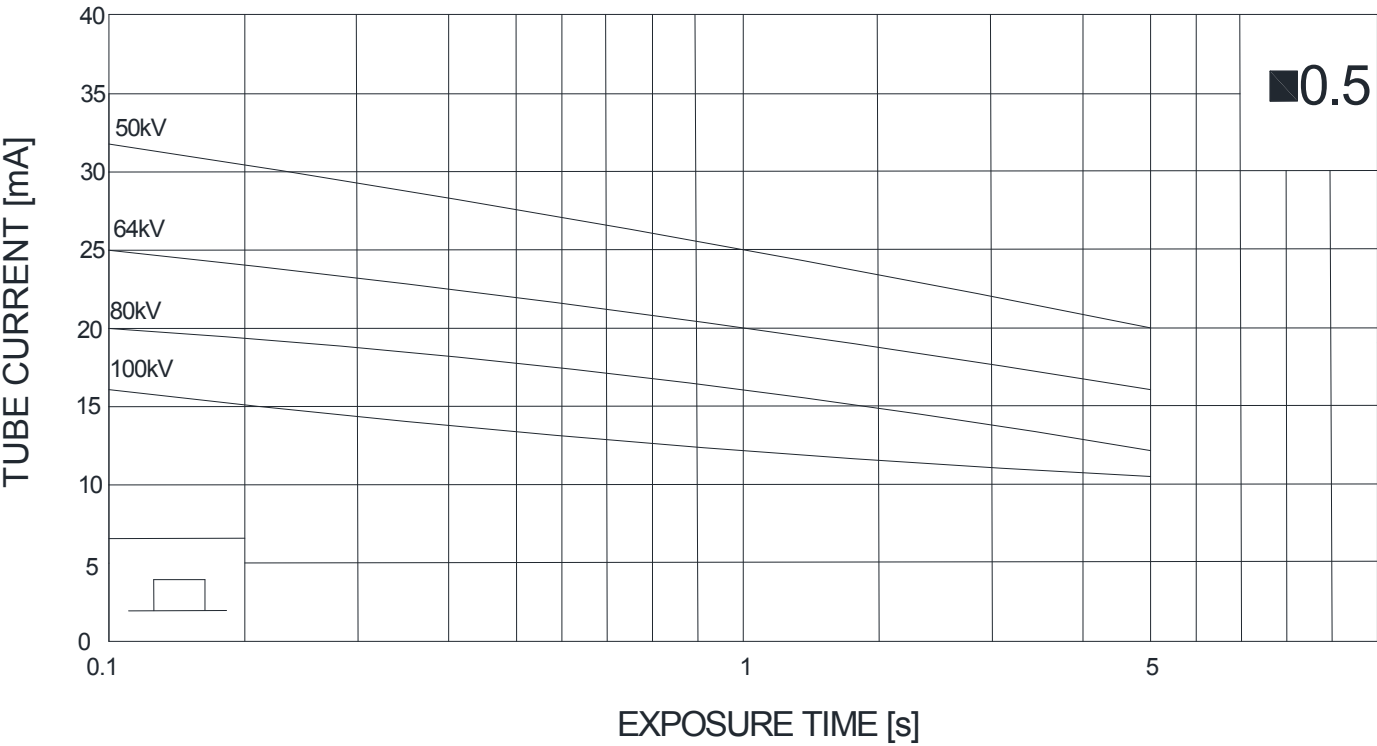

Emission & Filament Characteristics

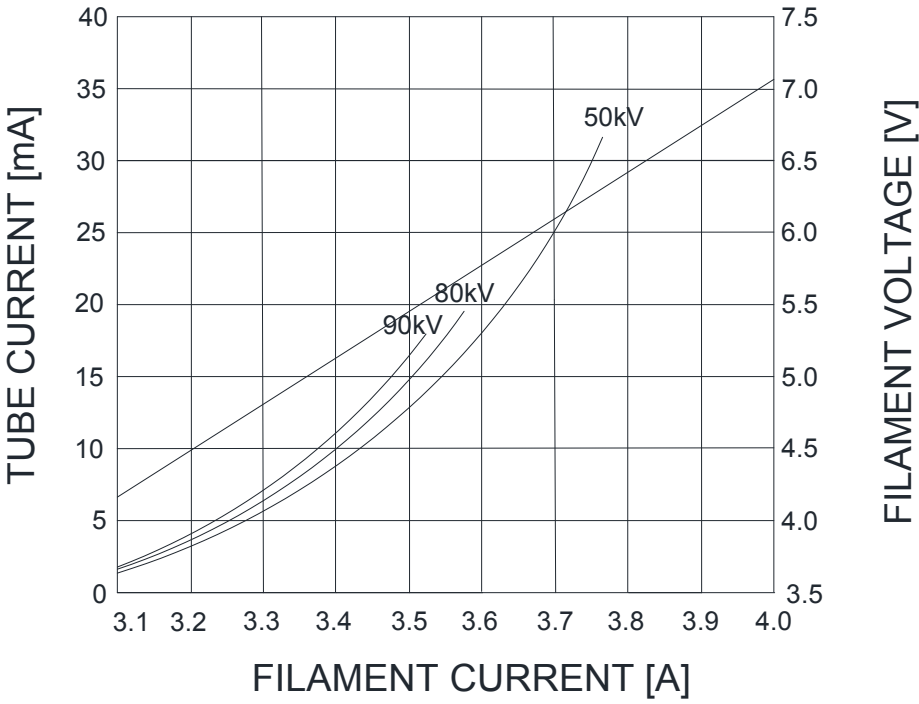


Anode Thermal Characteristics

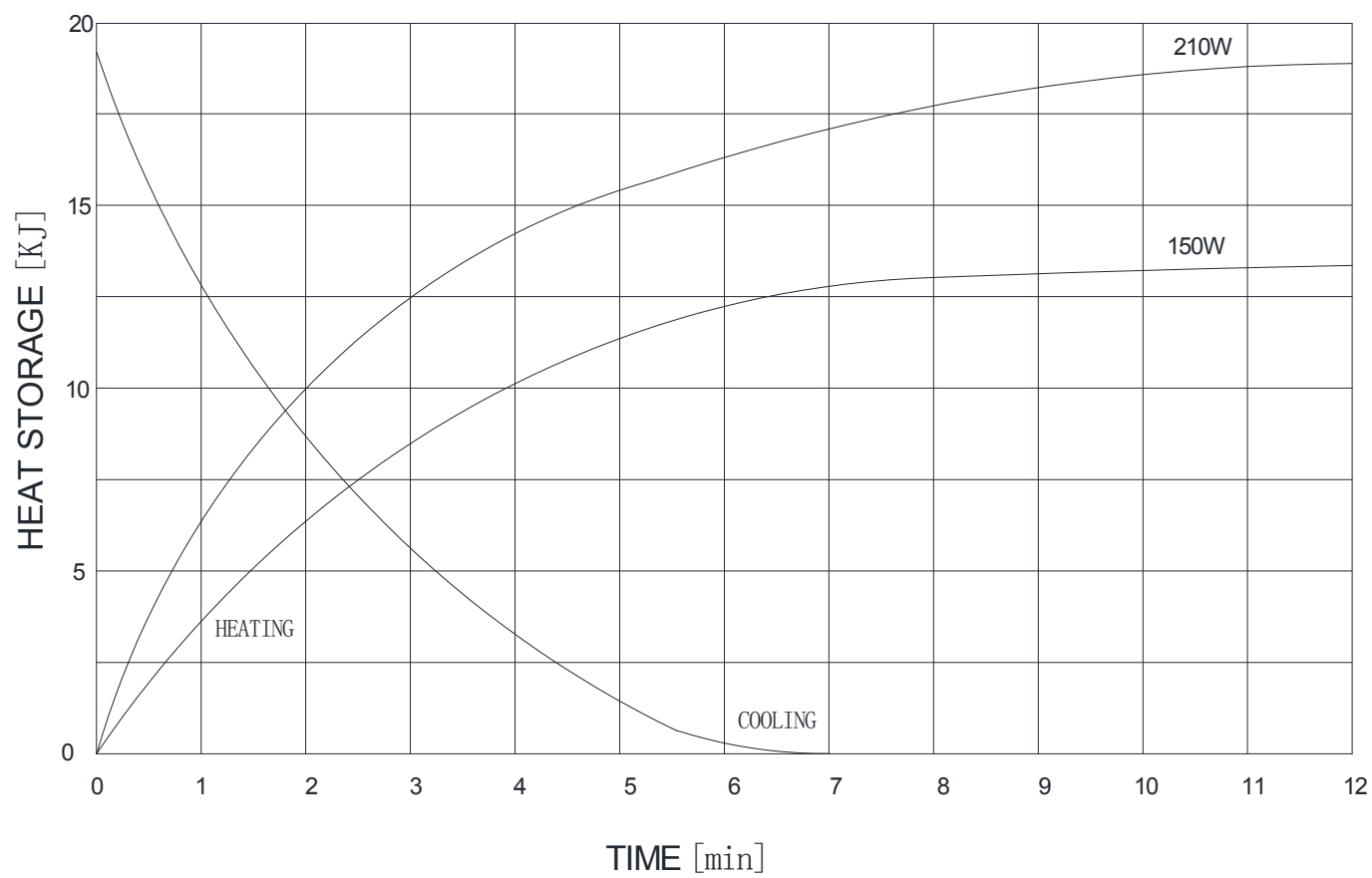


X-Ray Tube Dimensions -XD5-0.5-100 (12)

Unit: mm

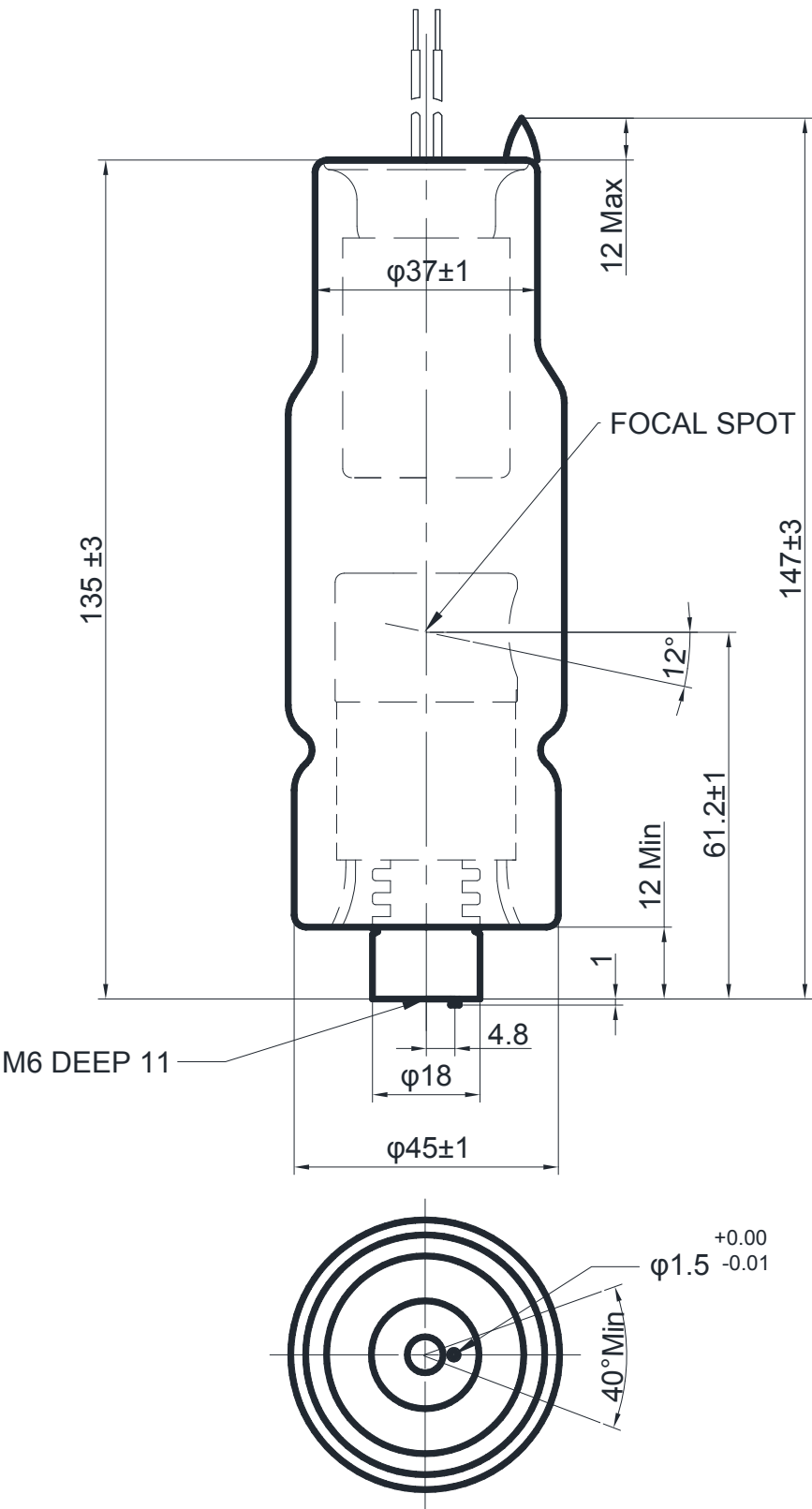


X-Ray Tube Dimensions -XD5-0.5-100S (12)

Unit: mm

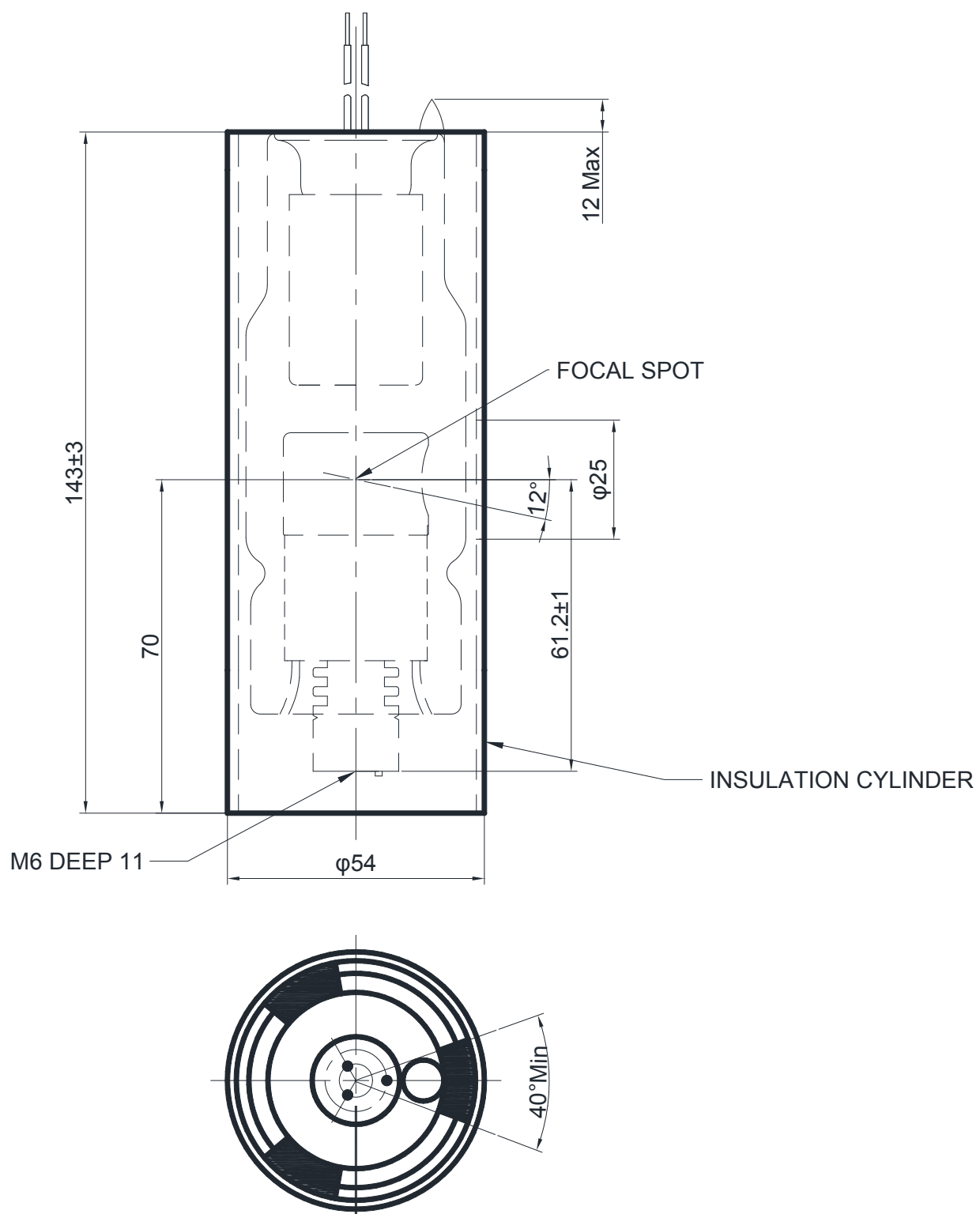

X-Ray Tube Dimensions -XD5-0.5-100SB (12)

Unit: mm

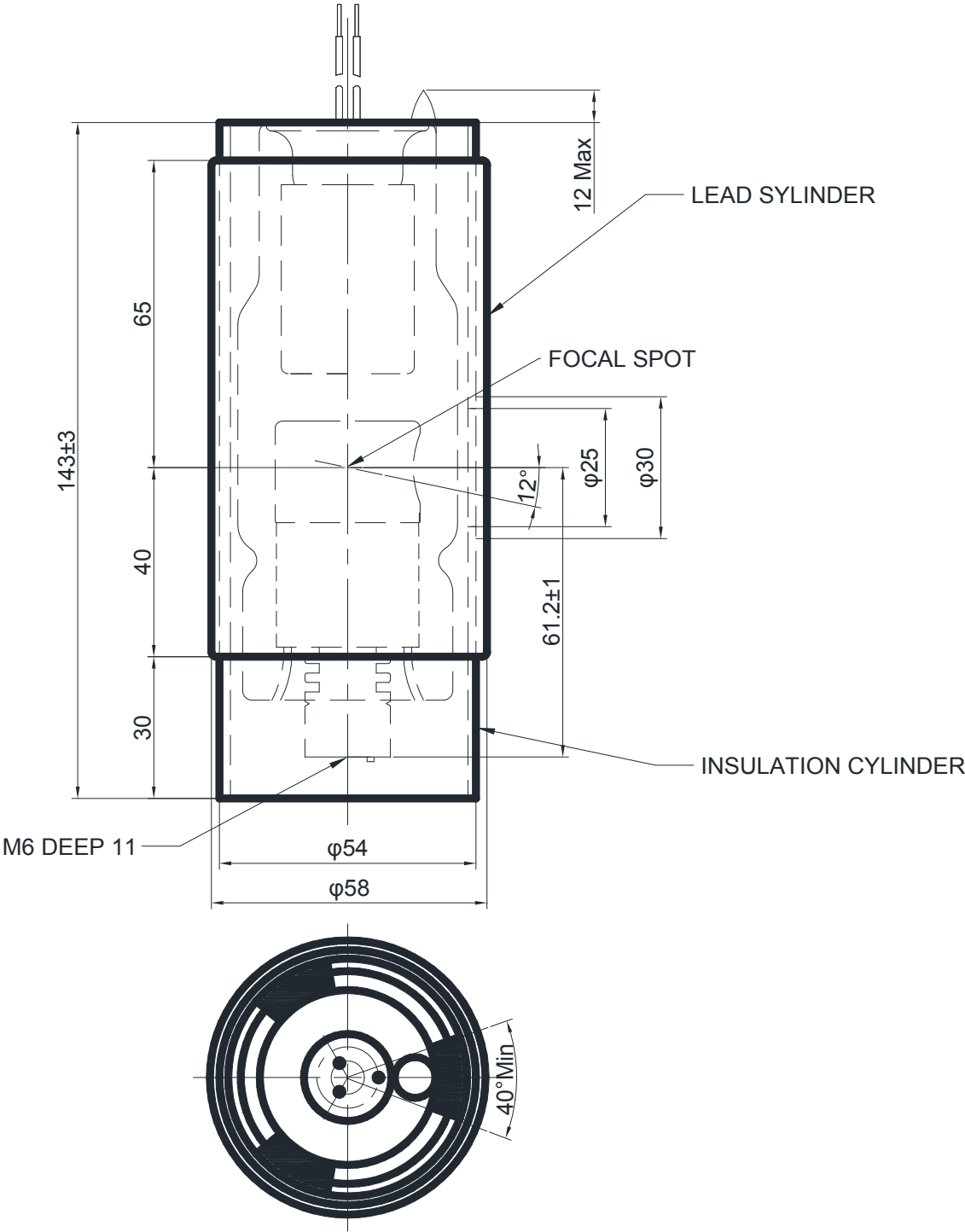

X-Ray Tube Dimensions -XD5-0.5-100SBR (12)

Unit: mm

